


# ENDOLINE® NASO-EXTEND

22

FULL EXTENSION

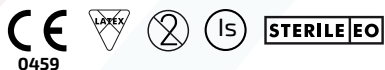


The **ENDOLINE® NASO-EXTEND** full extension drains bile or pancreatic secretions into an external receiving bag, and can also be used to rinse the nasobiliary catheter. This single use medical device is part of the **ENDOLINE® NASO-KIT** nasobiliary kit (for more information on the kit, please consult the data sheet page 18).

It has a 3-ways stopcock:

- 1 way to the patient, which is connected to the proximal end of the nasobiliary catheter via the Luer-Lock connector.
- 1 way for rinsing of the nasobiliary catheter, when required.
- 1 way to collect bile in an external receiving bag.

## AVAILABLE REFERENCE(S)

Reference	Length (cm)
AS0210013	135

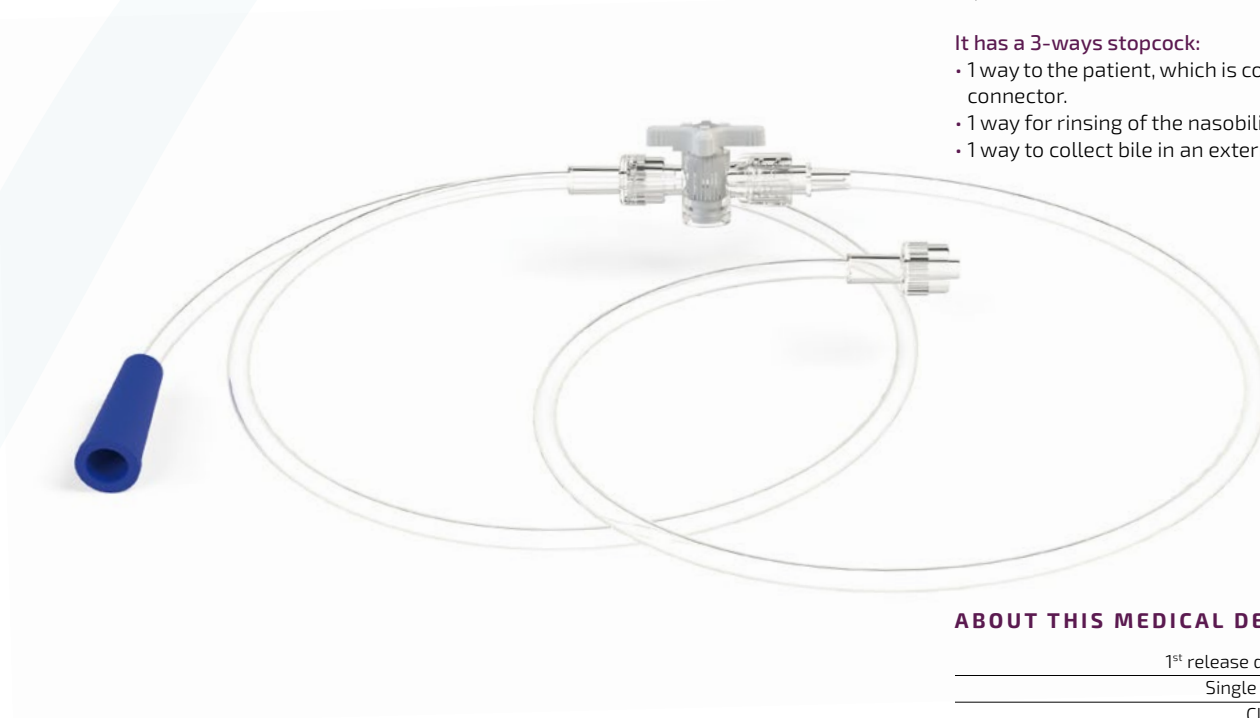

## ABOUT THIS MEDICAL DEVICE

1 <sup>st</sup> release date	1997
Single use	Yes
Class	Is
Sterile	Yes
Sterilisation method	Ethylene oxide sterilisation
Contains latex	No
Lifetime	5 years
Certification	CE
Packaging	Box of 5 units in individual sterile wrapping
Storage	In its original box in a cool and dry place

PRINCE MEDICAL SAS  
 ZA La Sente Du Moulin  
 64 rue Gutenberg - 60530 Ercuis - France  
[www.prince-medical.com](http://www.prince-medical.com)  
 Tel: +33 (0)3 44 26 54 76  
[prince-medical@omerin.com](mailto:prince-medical@omerin.com)

PRINCE MEDICAL  
 DEVICES

The information provided on this data sheet is indicative and subject to change without notice. PRINCE MEDICAL SAS shall not be held responsible for any incidents resulting from inappropriate use, failure to comply with the rules and regulations in force, or failure to follow the instructions for use provided. Medical devices sold by PRINCE MEDICAL SAS should only be used by a duly qualified practitioner or authorised person.

® Registered trademark of PRINCE MEDICAL SAS. Drawings and photos are not contractually binding. Reproduction is prohibited without the prior consent of PRINCE MEDICAL SAS.